

SEED STORAGE AND GERMINATION CHAMBER

Seed storage Chamber

Meditech is the largest manufacturer and exporter of seed storage chamber in India. Meditech seed germination chambers is used for seed germination, seed development and seed storage. These are special chambers which can maintain temperature, humidity, lights and CO2 levels.

Features of Seed Storage Chamber

- (1) These are designed exclusively for seed germination and seed storage and other research applications.
- (2) parameters such as Temperature, Light, humidity (optional) and CO2 (optional) can be controlled.
- (3) Touch Screen Controlled LED Controller with 12 types of alarms, RTC, Date and Calendar and data logging options
- (4) Program can be set in real and elapsed time.
- (5) Ramping mode is available for experiments
- (6) Multi profiles can be linked to create complex programs
- (7) SAIFIN program allows Seed Storage Chamber to access remotely.
- (8) Fluorescent Lamps are mounted on each Shelf vertically.
- (9) Light intensity upto 150 $\mu\text{moles/m}^2/\text{s}$ can be achieved of light irradiance measured @ 6" from lamps of seed storage chamber.
- (10) Programming is done using 7 Inch Touch screen controller

Airflow and circulation

- (1) Air is circulated horizontally throughout the chamber to maintain the desired temperature, humidity and CO2 levels.

Cabinet construction

- (1) Interior is constructed using 22-gauge 304 stainless steel.
- (2) Exterior is constructed using 18-gauge bacteria resistant powder coated CRCA aircraft steel
- (3) The Seed germination Chamber is completely self contained
- (4) The Thickness of cabinet is 3" (75 mm)
- (5) Two 2" diameter access ports are provided on Right side wall of the cabinet.
- (6) The Chamber floor is designed for drain which eventually end in a hose assembly.
- (7) The Chamber is supported with 4 castor wheels on with height adjustment facility.

Insulation

Insulation

- (1) The seed storage chamber is insulated using high density CFC free PUFF insulation.

Door

- (1) A single door provides full access to the chamber interior.
- (2) Gaskets with Magnetic strips provides a tight seal to door.

shelving

- (1) Fifteen tiers of Stainless Steel SS304 wire shelving

Refrigeration

- (1) The seed storage chamber is fitted with air-cooled refrigeration unit.
- (2) The refrigeration system allows easy access for maintenance and cleaning

Temperature range

- (1) 10°C to 44°C ($\pm 0.1^\circ\text{C}$) lights on and 2°C to 44°C ($\pm 0.1^\circ\text{C}$) lights off
- (2) Temperature controls, audible alarms, and visual indicators provided
- (3) Safety emergency Controls will shut down all the power to the chamber, and activates the alarms.

humidity control (optional)

- (1) humidity system for higher than ambient to 85% ($\pm 3\%$) lights on for set temperatures between 15° to 60°C
- (2) Humidity system for higher than ambient to 90% ($\pm 10\%$) lights off for set temperatures between 15° to 30°C

SEED STORAGE AND GERMINATION CHAMBER



The MTBBR Advantage

The Chamber series of equipment has been developed to provide the researcher and producers with dedicated technology for the development, storage and Germination of Seeds. Multi tiers of growing space maximize the chamber's footprint, making the MTBBR extremely space-efficient and ideally suitable where laboratory space is limited. Engineered and tapered air-handling plenums ensure consistent air-flow into the chamber. This, in combination with lamp configurations and canopies that are designed to maximize light dispersion, ensure a high degree of uniformity of all conditions within the growing area. This degree of uniformity has proven to be necessary for the efficient growth of modified forms of this species. All controlled parameters are managed and logged through Meditech's reliable 7" Touch screen controller giving the researcher maximum flexibility and documented results. Sound manufacturing principals, quality components and our ISO9001 quality control program ensures that this equipment performs as specified.



Lighting

The standard light intensity is 150 micromoles/m²/s using T8 fluorescent lamps or LED Lamps adjustable from 10% to 100%. Dimming ballasts provide precise control over light intensity. Meditech's own Closed-Loop Light Control maintains precise intensity of programmed levels.

Model	370 Litres	1000 Litres
Growth Range	3 growth areas (>1.1m ²)	8 growth areas (>3.3m ²)
Temperature Range	+10°C - +45°C	+10°C - +45°C
Temperature Fluctuation	±0.5°C @ 25°C	±0.5°C @ 25°C
Temp. Control Accuracy	±0.1°C	±0.1°C
Humidity Range	50% - 95% RH	50% - 95% RH
Humidity Uniformity	±5%	±5%
Humidity Control Accuracy	±1%	±1%
Humidifier Type	Ultrasonic Humidifier	Ultrasonic Humidifier
CO2 Control Range (optional)	Ambient - 5000ppm	Ambient - 5000ppm
Light Type	LED or FL T8 (blue, green, red, white)	LED or FL T8 (blue, green, red, white)
Light Intensity	up to 150 umols-1m-2 @15cm	up to 150 umols-1m-2 @15cm
Light Control	ON / OFF or 0% to 100% programmable	ON / OFF or 0% to 100% programmable
Program	Multi programs - multi steps	Multi programs - multi steps
Alarm	High and Low Temperature, Humidity alarm	High and Low Temperature, Humidity alarm
Data	Data Export by USB port	Data Export by USB port
Remote Control	Direct connection to unit by LAN to control HMI (by Internet and Wifi in option)	Direct connection to unit by LAN to control HMI (by Internet and Wifi in option)
Cooling System	Air cooled condensing unit / CFC free refrigerant	Air cooled condensing unit / CFC free refrigerant
Overview Dimensions	W890 x D720 x H1950mm	W1500 x D900 x H2100mm
Working Volume	370L	>1000L
Interior Material	Stainless Steel #304	Stainless Steel #304
Exterior Material	Coated steel / Stainless Steel	Coated steel / Stainless Steel
Insulation	PUF insulation	PUF insulation
Door	Insulated door	Insulated doors
Growth Area	3 growth areas (>1.1m ²)	8 growth areas (>3.3m ²)
Growth Height	25-30cm adjustable (2.5cm steps)	25-30cm adjustable (2.5cm steps)
Safety Device	Overheat protection and low temperature protection / Door open alarm / Over load protection / Compressor delay start protection	Overheat protection and low temperature protection / Door open alarm / Over load protection / Compressor delay start protection
Power	Adapted to the user electric system	Adapted to the user electric system

SEED STORAGE AND GERMINATION CHAMBER

Equipment meets electrical safety specification such as that of IEC (ClassI)

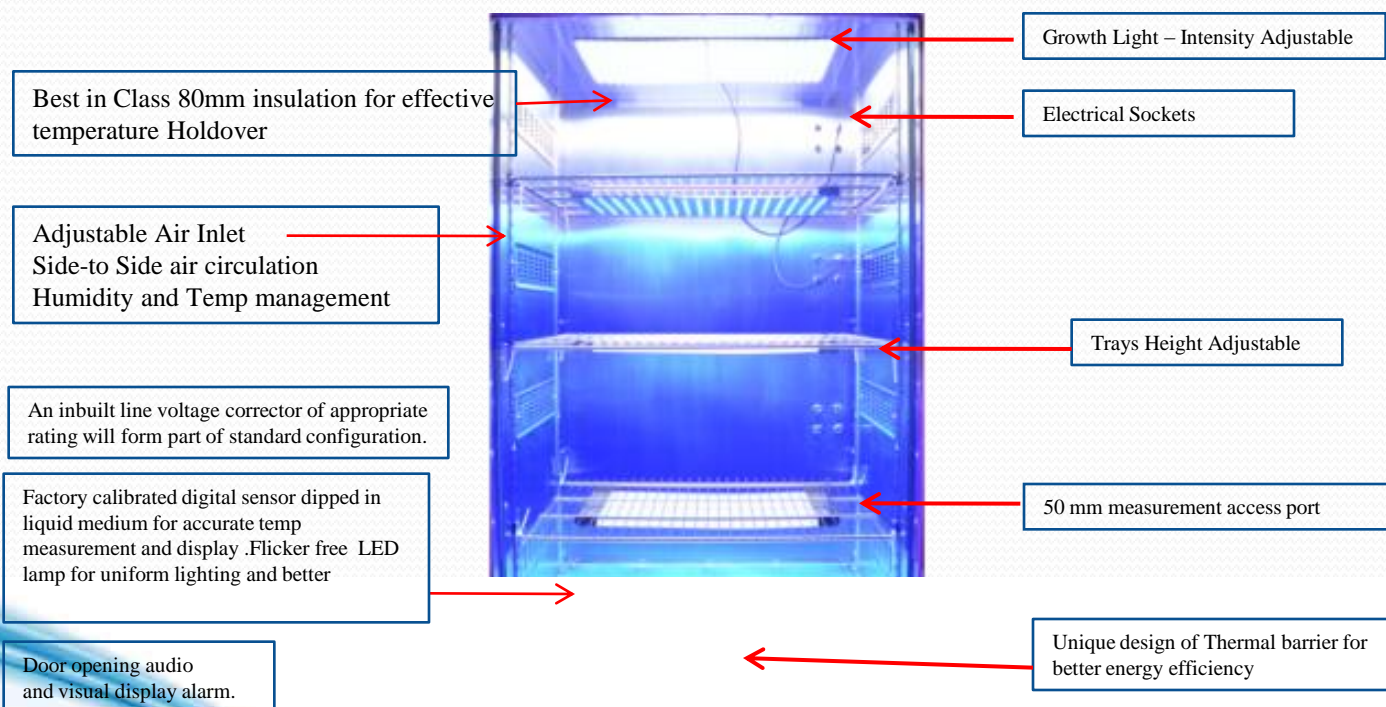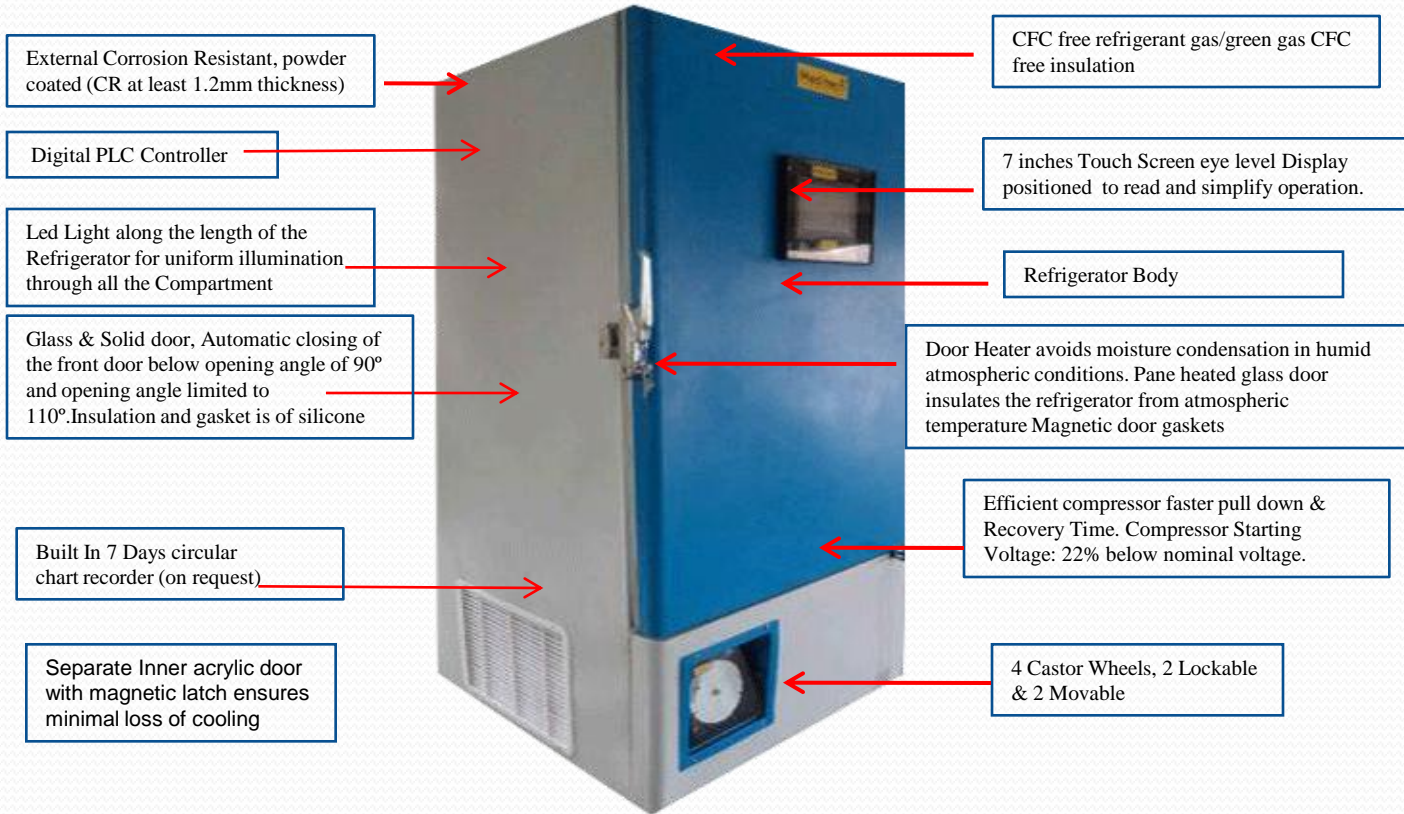

Technical Data

Standard Model(S): Inside S.S.304 with mirror polish & outside mild steel powder coated.

GMP Model(G): Inside and outside with SS 304 matt buff.



			Dimensions in meters					
Models	Temp	Capacity	Internal			External		
MTBBR	°C	ltrs.	W	D	H	W	D	H
MTBBR80	+10°C to +45°C	8000	2	2	2	2.16	2.95	2.16
MTBBR125	+10°C to +45°C	12500	2.5	2.5	2	2.7	3.45	2.16
MTBBR187	+10°C to +45°C	18750	2.5	3	2.5	2.7	3.95	2.7
MTBBR225	+10°C to +45°C	22500	3	3	2.5	3.16	3.95	2.7

Note: Other sizes are available on request.

SEED STORAGE AND GERMINATION CHAMBER

Web Based Software (SAIFIN53) - Remote Chamber and Room Management

Access Anywhere Anytime!! Around the Globe

Instantly connect and control your Chambers and Rooms without installing software

For more than a decade, Meditech SAIFIN53 software has gained popularity among operators and end users for its reliability and ease. SAIFIN53 is available remotely via a web browser using internet facility. SAIFIN53 is designed to meet the increasing demands for improved efficiencies. This web based software also adds extra capabilities of data logging, user-configuration and features of all alarm notifications via SMS and Email.

Features

- Eliminate the need for multiple software installations and Qualifications
- Web based software accessible through browser
- Compatible with server-based OS, android, Windows-7/8/9/10
- Ensure one-time software installation on your plants server and access the software anytime, anywhere through a web browser from any Desktop/Laptop (same LAN Network)
- View Status, take printouts or execute any of your tasks from anywhere
- Department wise segregation of users, equipment & Data
- Harmonization of reporting format (header, footer, etc) and other such operations through the plant

Centralized user defined Alarms, alerts and Notifications through SMS & emails

Configure the type & no. of alarms to be sent via email/SMS & configure delay time too

Running time calculation of MKT and Moving Average Component running hours to determine the shelf-life of major components

Real-time display of all connected equipment in single screen

Automatic data backup facility

System Diagnosis for troubleshooting via Team Viewer

Ethernet communication with chambers

Complying US FDA 21 CFR Part 11

Key Features

- A Truly innovative data handling and communication options
- Plug and Play fashion via Ethernet to connect local or plant networked PC or Laptop
- One can view LIVE system parameters and historical graphs as well as receive alarms, reports, data log files over email
- Settings with just a standard browser
- Access live or stored data remotely within the facility (LAN network) or from anywhere in the world



Secure audit trails
-Audit protection
-Provide certified proof of process control

FEATURES AND FUNCTIONS

A REVOLUTIONARY MONITORING & ALARM SYSTEM

Say No To Manual Recording

- Error Free
- Separate signatures for data reviews & data approval
- No mistake in record maintaining
- Data analysis
- Tabular and graphical data representation



Wireless Display with Humidity, Temperature CO2 Sensors for Visibility and audit

Record and alarm Temperature, Relative Humidity, Flow, Differential Pressure, CO2, O2, LN2, Status and Contacts

USES : Freezers, Stability Chambers, Warehouse, Labs, Cold Rooms and Blood Banks

View Multiple areas or selected locations on PC



View Data from Any part of the World
-Browser base application software
-view from PC, smart phones and tablets



Alarms
Alarms for parameter deviation
SMS Alerts
Alarm log with acknowledgement and user comments
Email alerts for supervisor


Audio and Visual Alarms for Settable Range

Compact size chart recorders and alarms system

Data recorders are easy to relocate and swap for service and calibration



Receive SMS and Email alarm notification

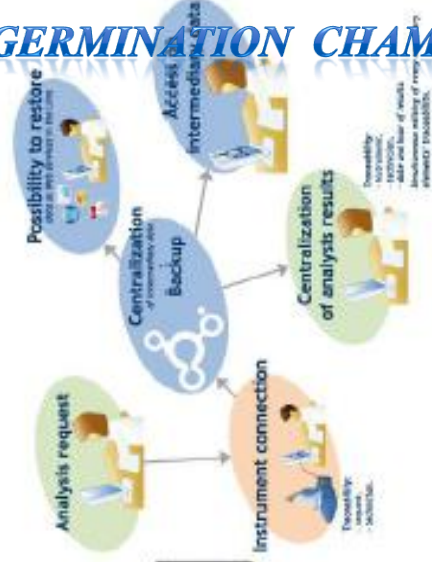

Traceability

- Instant response to data request
- Continues Data saving
- Data Protection and backups



Supervisor Function

Alerts users for pending activity
Escalates unattended issues to higher authority

SEED STORAGE AND GERMINATION CHAMBER



Screen Shots of Touch Screen Controller










ISO 9001:2008, ISO 14001, ISO 13485
WHO-GMP & CE Certified Company

Meditech Technologies India Pvt Ltd


No.6 President Terrace.Takli Road, Dwarka
Nashik – 422011, India
Tel No: 0253-6418632
info@newmeditech.com;huzaifa@newmeditech.com

 16/02/15		PARAMETERS		 17:50:31	
DESCRIPTION	VALUE	UNIT	DESCRIPTION	VALUE	UNIT
Heater On Time	±99999	Min			
Temp Set RTD1	±9999.9	°C			
Temp Set RTD2	±9999.9	°C			
Comp. Sel.	Comp1				
Comp 1 Time	±9999.9	Min			
Comp 2 Time	±9999.9	Min			





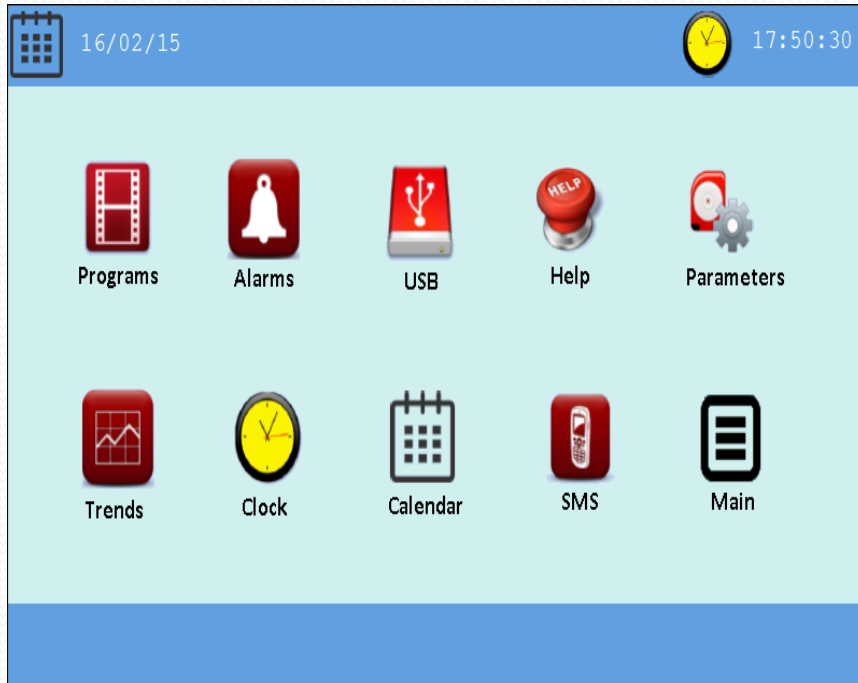

 MENU

SEED STORAGE AND GERMINATION CHAMBER


Screen Shots of Touch Screen Controller




SEED STORAGE AND GERMINATION CHAMBER

Screen Shots of Touch Screen Controller

The screenshot shows the SMS application interface. At the top, there is a status bar with a calendar icon, the date '16/02/15', the text 'SMS', a clock icon, and the time '17:50:32'. The main area has a light blue background. It contains a red label 'ENTER MOBILE NO:-' followed by a text input field with a black border and a black mask '#####'. Below this is another red label 'SEND SMS:-' followed by a red button with a white mobile phone icon. A horizontal line separates this from a section with a white button labeled 'RECEIVE MESSAGES'. Below the line is a large text input field with a blue border and a blue mask '#####'. At the bottom, there is a dock with six icons: a red square with a white graph, a red square with a white film strip, a red square with a white bell, a red square with a white mobile phone, a red circular button with a white 'HELP' label, and a red square with a white USB icon. On the far right, there is a black square icon with three white horizontal lines and the text 'MENU'.








 16/02/15

ALARMS - REALTIME

 17:50:33

ALARMS	ON DATE TIME	OFF DATE TIME	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N
XXXXXXXXXX/YY HH:MM:SS	DD/MM/YY HH:MM:SS	DD/MM/YY HH:MM:SS	Y/N

ACK
ACK ALL
HISTORICAL








MENU