


Rabbit Restrainer (Stainless Steel)

Meditech Technologies India Pvt. Ltd is the largest manufacturer and exporter of Rabbit Restrainer (Stainless Steel).

These restrainers are made of stainless steel 304. The adjustable neck portion is made of polished Bakelite or teakwood material.

The length can be adjusted by moving the S.S. plates and then tightening by using wing nuts. The neck can be fixed by opening the top S.S. plate by releasing the lock and then adjusting the height of upper wooden block and then again locking the S.S. clamp.

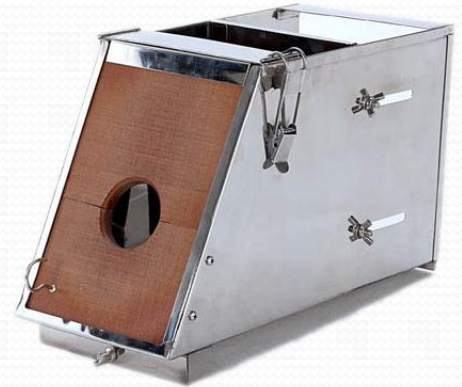

Ordering Information

Model No.	Animal Weight Range	Overall Dimensions (Mm)	Restraining Min Height (Mm) From Bottom	Restraining Max Height (Mm) From Bottom
Lar1s	3 To 7.5 Kg	380 X 175 X 230	220	285

Fully Adjustable Acrylic Rabbit Restrainer

Meditech Technologies India Pvt. Ltd is the largest manufacturer and exporter of Rabbit Restrainer (Stainless Steel).

These restrainers are made of stainless steel 304. The adjustable neck portion is made of polished Bakelite or teakwood material.

The length can be adjusted by moving the S.S. plates and then tightening by using wing nuts. The neck can be fixed by opening the top S.S. plate by releasing the lock and then adjusting the height of upper wooden block and then again locking the S.S. clamp.

Features

- Incorporates 2 designs in single instrument
- Ease of use and access to the rodent
- Made from thick perspex material



Ordering Information

Model No.	Animal Weight Range	Overall Dimensions (mm)	Design Usage Style	Restraining Min Height (mm) from bottom	Restraining Max Height (mm) from bottom
LARIA	3 to 7.5 Kg	380 x 175 x 230	Style I (curved back plate)	220	470
			Style II (2 plate system- Curved top and straight back plate)	210	470